



CUBE 30 ATV

A NEW APPROACH TO NON-INVASIVE VENTILATION

TECHNICAL SPECIFICATIONS

Ventilation modes

CPAP
PSV S / BiLevel S
PCV / BiLevel T
PSV / BiLevel ST
APCV

Pressure range

IPAP 5-30 / EPAP 3-25
cmH2O CPAP 3-20 cmH2O

Pressure support

0-25 cmH2O

Target Volume:

OFF, 200 - 2.000ml

Target Minute Volume:

OFF, 2.000 - 15.000ml

Respiratory rate

4-30 bpm

Inspiration time

0.3 - 5.0 seconds

Inspiration Trigger

Five Level settings

Expiration Trigger

Auto, Manual settings 70%-30%

Trigger Lock

0.3 - 8.0 seconds

Alarms

11 System/Therapy alarms

Power supply

Internal Switch Mode
Power Supply
100-240 V, 50/60 Hz

Dimensions

Size (W x H x D)
168 x 108 x 182mm
Weight 1,6 Kg


**Made in
Germany**

Maximum Performance, minimum footprint & Made in Germany!

The NEW Cube 30 ATV from Drive DeVilbiss International is an innovative non-invasive volume supporting ventilator. The Cube 30 ATV incorporates intuitive state of the art functionality allowing safe and easy application of the device in a homecare or clinical environment. With its unique trigger technology and target minute volume algorithm, the CUBE 30 ATV will optimise therapy to meet individual patient requirements.

DESIGN

- The device includes a large 4.3" coloured screen display incorporating easy to use tile design menu screens for clear convenient user navigation.


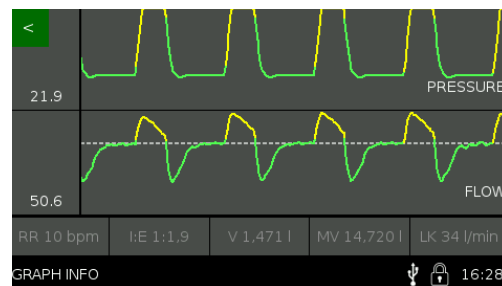


TECHNOLOGY

- High sensitive flow based trigger detection with trigger failure auto-detect function.
- Contemporary volume assured pressure support which can be set based on breath-by-breath target volume or target minute volume.
- Extensive device alarm management, including on screen display alarm protocol.
- Swiss made turbine blower promising high quality performance for extended, reliable operation.


INFORMATIVE

- Showing clear simple text information the device displays corresponding data settings including user notifications for incorrect parameter settings.
- Three different menu screen colours identifying Patient and Clinical menu modes.

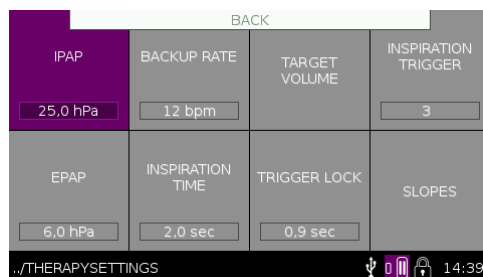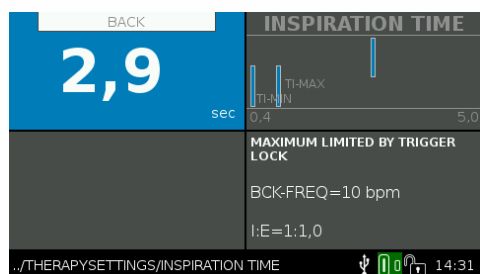

VARIOUS ADJUSTMENT OPTIONS

- Two independent adjustable ventilation settings allowing patients to change therapy mode for daytime and night time use.
- The Mask Check function ensures a correct fitting mask even before therapy starts. On screen values identify leakages in order for the mask to be adjusted and corrected.

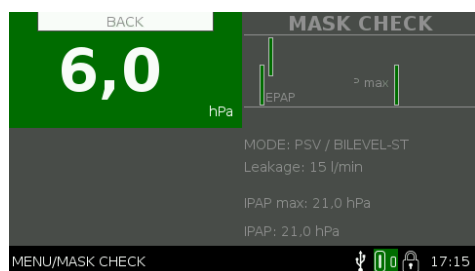

OPTIONAL HUMIDIFIER

- User friendly and high humidity performance for maximum comfort during therapy.



- The integrated humidifier connects underneath the Cube 30 ATV without increasing the foot print of the device and therefore saving valuable space.