



HbO-3000

There's a new kid on the block!

SPECIFICATION

Display Type	OLED
Dimensions	57 mm x 31.5 mm x 30.5 mm
SpO2	Range: 35% - 100%. Accuracy: 70% - 99% \pm 3%
Pulse rate	Range: 30 – 250. BPM Accuracy: 30 – 250 BPM \pm 2%
PI Range	0 – 30% 0.1% Accuracy
Power	2 x 1.5 AAA batteries (20 hours use)



Drive DeVilbiss Healthcare would like to introduce the next generation device in the HbO Pulse Oximeter range.

The HbO-3000 fingertip pulse oximeter contains all the necessary features you would expect for SpO2 spot checking ... and more!

- Integrated SpO2 and pulse rate alarms notifying the user of high and low values
- Perfusion Indexing (PI) identifying the best location for signal detection
- 5 different brightness settings and 6 screen display option for user preference
- Low power requirements – 2 x AAA batteries
- Auto off function to reserve battery life
- Designed and built to withstand the challenges of daily use!

PRODUCT FEATURES

The HbO-3000 features an optional six directional display mode for operator convenience, Perfusion Indicator, 5 OLED brightness levels, and user defined parameter settings for both SpO2 and pulse rate criteria indicating visual & audible alerts for high and low values.

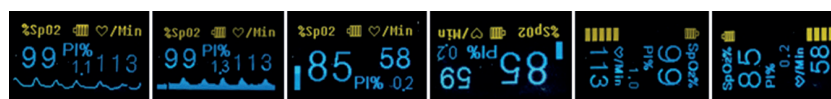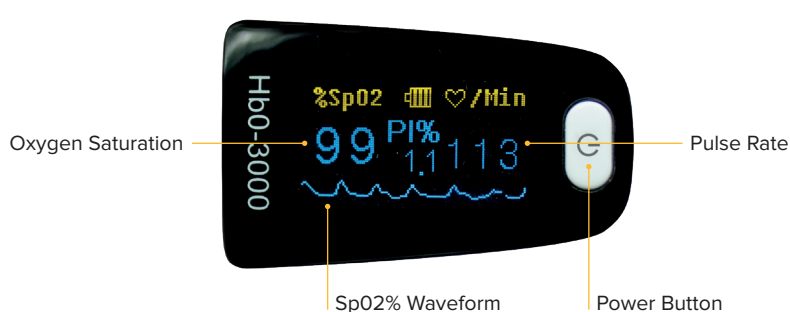

SpO2, Pulse rate and Pulse rate bar, Plethysmograph waveform, Perfusion Index and battery status indicator

Drive DeVilbiss International

DeVilbiss Healthcare GmbH
Kamenzer Straße 3
68309 Mannheim, Germany

UK +44 (0) 121 / 521 31 40
DE +49 (0) 621 / 1 78 98 0
kontakt@devilbisshc.com

www.drivedevilbiss-int.com